

# [ExerCube] — Recorded Metrics, Including Descriptions

Markdown conversion of `ExerCube Data.docx` (Sphery metrics documentation, received July 2026). Definitions below are Sphery's own wording.

## SpeedCage program

Performance metrics available for both the entire workout and individual laps.


*Tablet:*



*App:*



- **Score:** Total points scored based on targets hit
- **Distance:** Distance covered in metres
- **Calories burned:** Number of calories burned
- **HR**
  - **HR Average:** Average heart rate
  - **HR-Max:** Maximum heart rate achieved
  - **HR-Recovery:** Difference between HR value at the start and end of a break
  - **HR-Zones:** Time spent in zones calculated for the individual (image of cube wall; above/in app for more detail)

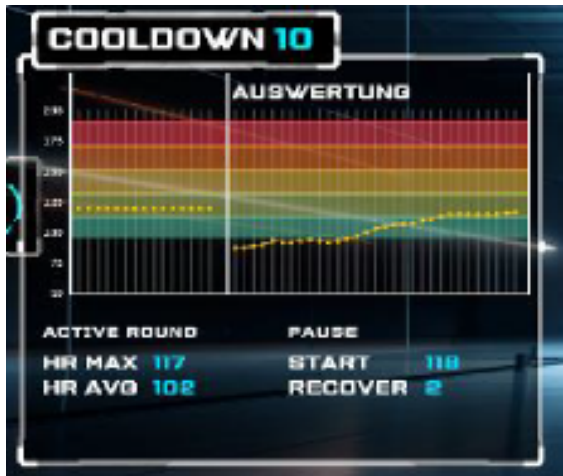

- **Reaction Time**


- **RT-Average:** Average reaction time (time from appearance to hitting the target)
- **RT Analysis:** Distribution of reaction times across thirds of the target display time (0–4 s)

- **Targets**

- **Analysis:** Number of normal, bonus, and incorrect targets hit
- **Distribution:** Analysis of targets touched per wall (totals 100% if all targets are hit)

- **Heat maps**

- **Targets (hits):** Detailed analysis of correct targets hit per field (left image)
- **Targets (reaction):** Detailed analysis of reaction times for correct targets hit per field (right image)

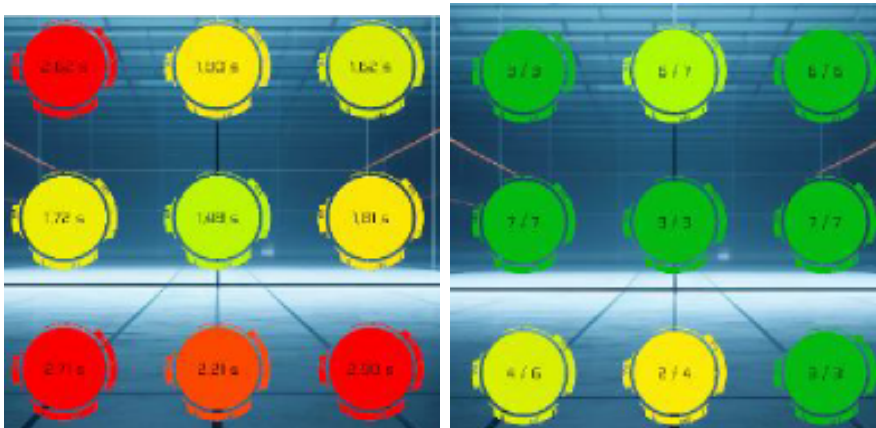

- **Brain Speed Assessment** (*available once a month*) Tracks improvements in reaction time using a standardized training protocol with challenging cognitive tasks and a fixed benchmark round for direct performance comparison.



# Racer program

Tablet:



App:



- **Points:** Score achieved
- **Body Score:** Percentage of exercises performed correctly (%)
- **Brain Score:** Percentage of timings performed correctly (%)
- **Dual-Flow Score:** Combination of Body and Brain Scores (%)
- **Biggest Combo:** Highest number of exercises performed correctly in succession
- **Burned Calories:** Number of calories burned
- **HR-Max:** Maximum heart rate achieved
- **HR-Average:** Average heart rate
- **Distance:** Distance value converted to moderate jogging based on workout duration, number, and type of exercises

# Further information

## **Assessment program (*no reporting currently available*)**

- **Simple Reaction:** Reaction speed
- **N-Back:** Reaction, concentration & memory
- **Trail Making:** Cognitive flexibility, processing speed, attention, reaction and overview